


Trust M.A.S. Medical with all your Digital Full Body and Dental Imaging Needs



scanX QUANTUM

A Unique digital advantage for your practice. Full body and dental imaging in one.

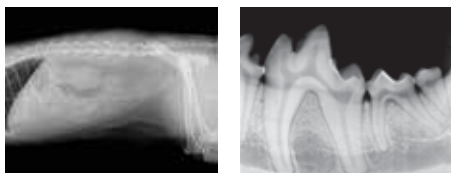

The ScanX Quantum is the only digital imaging system capable of producing diagnostic quality digital images in all standard sizes up to 14" x 17" and all dental sizes. ScanX Quantum is easy to learn and operate, and set-up is quick and simple. ScanX Quantum produces diagnostic quality digital images in just seconds. Imagine the time and efficiency you will add to each day with an upgrade to digital with the ScanX Quantum. ScanX Quantum is the perfect solution for veterinary practices looking to upgrade to digital and increase revenue by adding dental X-ray examinations to their routines.


With the purchase of any Digital x-ray system receive a FREE laptop valued at \$899.00



scanX duo

DIGITAL IMAGING SYSTEM

Our smallest ScanX at less than one cubic foot.

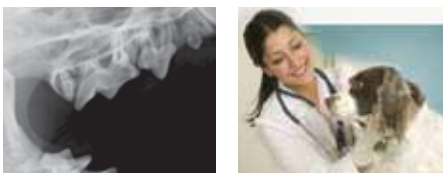

Now you can have the same great quality as our other ScanX products compacted to fit almost anywhere in your office.

The exclusive Independent Track Control (ITC) feature allows the processing slots to be simultaneously or independently controlled. A single ScanX Duo can be operated by two users at the same time with different patients making it perfect for use in a shared X-Ray pass-through.

ScanX Duo uses Phosphor Storage Plates that are ultra thin, wireless and flexible allowing for accurate positioning within the oral cavity. Using PSP's minimizes image elongation, foreshortening and retakes. The Plates are very inexpensive.

As with the original ScanX, the Duo uses In-Line Erase to read the PSP and automatically erase it in the same continuous cycle; simplifying the image process.

Compact, Affordable, Versatile Digital Imaging.



FireCR

The CR system that makes your job easier.

Increase the productivity and professionalism of your veterinary practice with the first full-featured computed radiography system offering flexibility that lets you grow.

The next-generation FireCR+ Veterinary CR Reader performs equally well in busy veterinary clinics scanning images of household pets as it does when transported to a stable, ranch, or zoo for examinations of large animals.

A Practical CR Solution

The compact, lightweight FireCR+ Veterinary Reader can be wall-mounted. Need portability? Even when the unit is wheeled around a clinic or transported for on-site imaging of large animals, auto calibration technology ensures consistently high-quality x-ray images.

The CR reader is available in several models, which can be affordably upgraded in the field to meet changing needs:

- FireCR+ Veterinary – **60 CR Reader** for veterinary hospitals and busy clinics
- FireCR+ Veterinary – **80 CR Reader** for horse exams, which require numerous x-rays
- FireCR+ Veterinary – **40 CR Reader** for small-animal veterinary clinics
- FireCR+ Veterinary – **20 CR Reader** for lower-volume small-animal veterinary clinics

Your Number 1 Canadian Distributor for ScanX Digital Imaging Systems, Parts and Accessories

- **Phosphor Storage Plates for Dentistry**
- Size #0 / #1 / #2 / #3 / #4 and our new size #6 now available
- Barrier Envelopes
- Plate Guides
- Plate Transfer Box
- Cleaning Sheets
- Cleaning Wipes
- **Phosphor Storage Plates for Full Body**
- Plate / Plate Protector and Cassette 8 x 10
- Plate / Plate Protector and Cassette 10 x 12
- Plate / Plate Protector and Cassette 14 x 17
- Cleaning Sheets

IT Support and ScanX Software Upgrades

NOW Available your own Preventative Maintenance Service Agreement

If you already have a ScanX System over one (1) year old, please read this:

Like all precision products, the ScanX requires a certain amount of care. A well-organized maintenance program aids dependable equipment operation and reduces problems to a minimum. Routine checks help to detect general overall wear, and replacement of parts can often be made before problems occur.

Understanding this, we at M.A.S. Medical with Allpro Imaging have established two basic maintenance kits that will help ensure continuous operation of the ScanX Digital Imaging System. Please give us a call for more details.



From left to right: Phosphor Storage Plates, Barrier Envelopes, Plate Guides, Plate Transfer Box



36 Lauraview Crescent, Unionville, Ontario, L6C 2R4

www.masmedicalsupplies.com

Email: sales@masmedicalsupplies.com Phone: (905) 927-0548 Fax: (905) 887-5306