




GS MEDICAL PACKAGING INC.

30 Years of Advanced Solutions in all of your Infection Control Needs

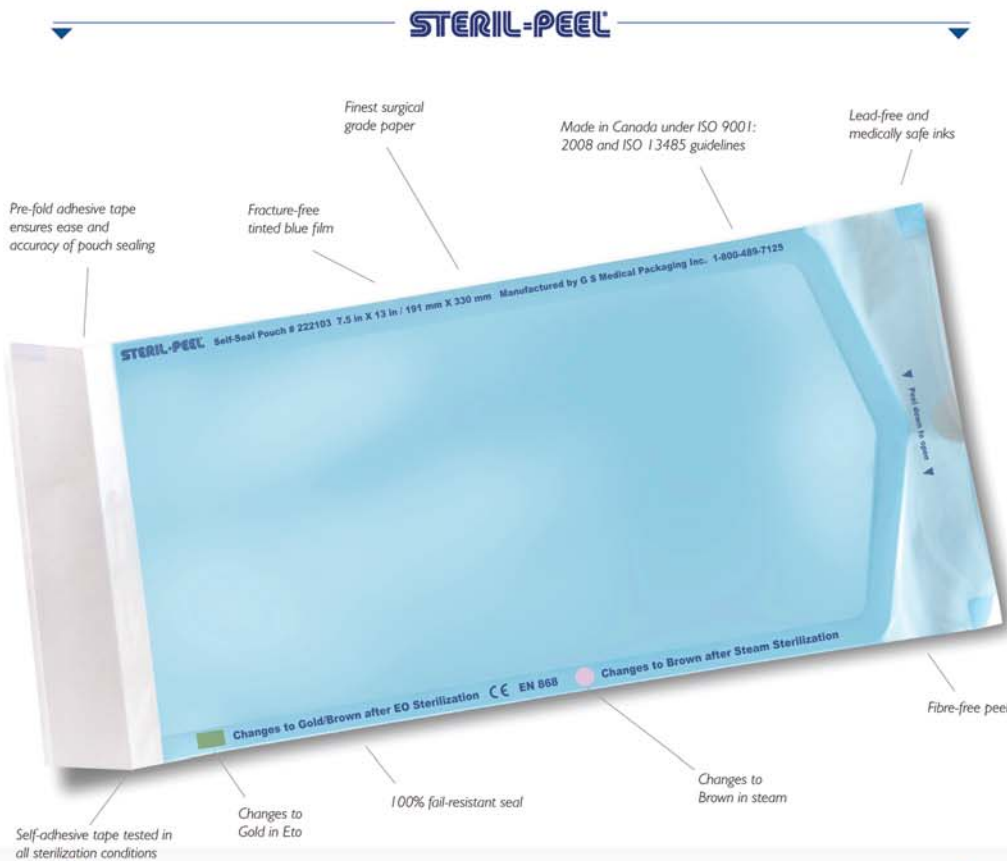

All STERIL-PEEL* product are manufactured from medical grade bleached kraft paper. Hydrogen peroxide, rather than chlorine is used as the bleaching agent.

They are made from virgin pulp, ensuring optimum consistency of the paper sheet, providing superior barrier properties.

STERIL-PEEL* self-sealing pouches are manufactured with a pre-fold (see diagram) to ensure ease and accuracy of pouch sealing.

The adhesive tape is made of a polyester film coated with medical grade acrylic adhesive, designed for use in steam and ethylene oxide sterilizers.

QUALITY MANUFACTURING

- Compatible with all of the most commonly used sterilization methods.
- Validated by independent laboratories to ensure that all parameters required have been met in compliance with ISO standards.
- Easy identification of contents through transparent blue film. Blue film vs white film – Blue film being a stronger identifier of defects, punctures and/or tears in film.

Water-based adhesives
Lead-Free Inks
Dual Indicators
Latex-free

GS Medical Packaging enjoys over 30 years of being the pioneer in high quality medical pouches and sterilization products. We are one of the few companies that manufacture self-seal pouches with Class 4 indicators. Since our inception, we have created packages that are a cut above the rest in strength, sterility and protection for medical devices, components, trays and cassettes.

Our superiority in design, quality and function is indisputable. These are trademarks of our brand.

GS Medical Packaging's reputation continues to grow both Nationally and Internationally due to our exceptional products and our extensive in-house knowledge of infection control.

Infection Control Above the Rest



Item No.	Size	Millimeters	Case Quantity
222100	3.5"x5.25" (89 mm x 133 mm	5 x 200
222101	3.5"x9"	89 mm x 229 mm	10 x 200
222102	5.25"x10"	133 mm x 260 mm	10 x 200
222103	7.5"x13"	191 mm x 330 mm	5 x 200
221104	5.25"x15"	133 mm x 381 mm	5 x 200
222816	8"x16"	203 mm x 406 mm	5 x 200
222555	12"x15"	305 mm x 381 mm	2 x 200
222565	3.5"x22"	89 mm x 559 mm	5 x 200
222218	12"x18"	305 mm x 457 mm	2 x 200

STERIL-PEEL



Professional's Choice*
Steril-Sure*



9001
13485
22000



PROUDLY DISTRIBUTED BY:



M.Q.S. MEDICAL SUPPLIES INC.
MEDICAL • DENTAL • VETERINARY

Contact us for more information:
www.masmedicalsupplies.com

Email: sales@masmedicalsupplies.com
Phone: (905) 927-0548 Fax: (905) 887-5306



GS MEDICAL PACKAGING INC.

PROUDLY CANADIAN



Your most trusted source
for infection Control Supplies